


CRS Electronic Crimpers and Decappers

Hand-held Rechargeable Crimpers and Decappers



- Precision, rechargeable crimpers and decappers for 8, 11, 13 & 20 mm caps.
- Easy to use; comfortable grip.
- Hundreds of reproducible crimps from one charge.
- Adjustable crimping force.
- Includes standard Black & Decker VersaPak® Gold battery and charger.

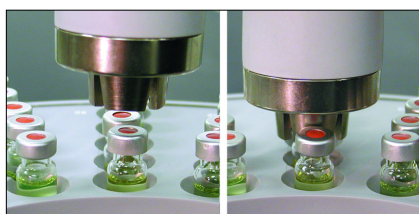

The Electronic Crimper is designed to fit around vials placed in standard autosampler trays. The adjustable setting provides a precision crimp, vial after vial.



11 mm Electronic Decapper



The cycle of these tools are controlled precisely with an internal counter.




Crimp headspace vials with the push of a button.



Adjustment buttons for precise control.

The CRS electronic crimping tools are powered by a standard Black and Decker Versapak® battery and will cycle hundreds of times on a single charge. The cycle of all of these tools is controlled precisely with an internal counter and is adjusted with two buttons on the side of the case. The tools fit comfortably in the hand and weigh about 600 grams. When used in a seated position the jaws can be positioned easily around closely-spaced vials in standard autosampler trays.

The complete tool comes with a Versapak® Gold rechargeable battery, which uses environmentally friendly NiMH technology, and a charger. Recharging generally takes 6-9 hours. The CRS electronic crimping tools are available in versions for most countries.

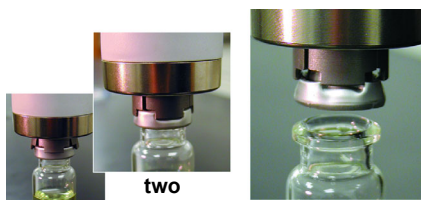

one

two

three!

Decapping 20 mm caps is now easy.

Description

Item No.

Electronic Crimper, 8 mm, Complete Kit*	
Electronic Crimper, 11 mm, Complete Kit*	
Electronic Crimper, 13 mm, Complete Kit*	
Electronic Crimper, 20 mm, Complete Kit*	
Electronic Decapper, 11 mm, Complete Kit*	
Electronic Decapper, 13 mm, Complete Kit*	
Electronic Decapper, 20 mm, Complete Kit*	

Accessories

Replacement VersaPak® Gold (NiMH) Battery

* Complete kit includes tool, one battery and battery charger.

